



Multichannel Cell Counter CZI Cellpose


Quick User Guide

Multichannel Cell Counter CZI cellpose automates cell detection and counting in multiple single plane files (no z-stacks). For each image, single cells are detected by cellpose (deep learning-based segmentation). For each channel and cellular mask, a staining is considered positive if a minimum number of pixels are above a given threshold. Combinatorial filters for cell counting can be defined based on staining (e.g., A+ B+ C-). Parameters are best extracted from FIJI pre-analysis.

Step 1 Detect Cells

1. Set Nuclei **Channel** (check ZEN for DAPI channel)
2. Set Nuclei **Diameter (pixels)** – use FIJI for nuclei area measurement in pixels
3. Choose **Cellpose model** (e.g. cyto2 for cells)
4. Set **A**, **B**, and **C** to Discard (D)


5. Select **Test file**
6. Press **GO!**

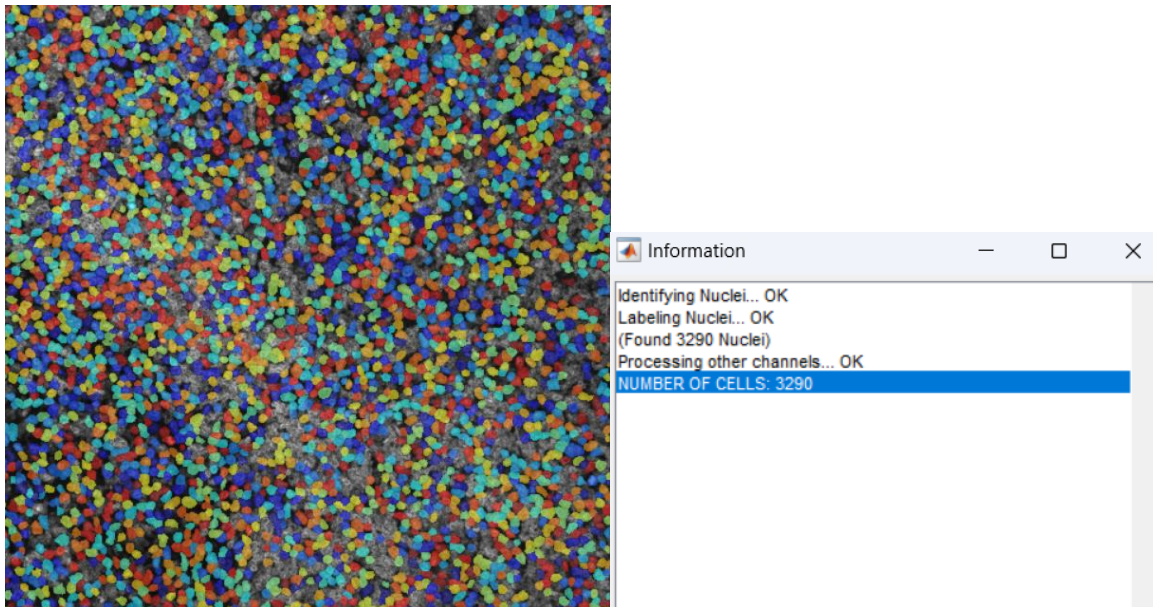


1 - Output image with nuclei detected by cellpose model cyto2

Step 2 Detect Positive Cells

1. Set A, B and C **channels** (check ZEN for correct numbering)
2. Set A, B and C **Threshold** values (use FIJI for threshold adjustment)
3. Set **Minimum Pixel Number** – the minimum number of pixels above threshold for positive staining
4. Set A, B, and C to either Y (Yes), N (No) or D (Discard) – set Discard for non-existing channels
5. Press **GO!** to detect positive cells in test file

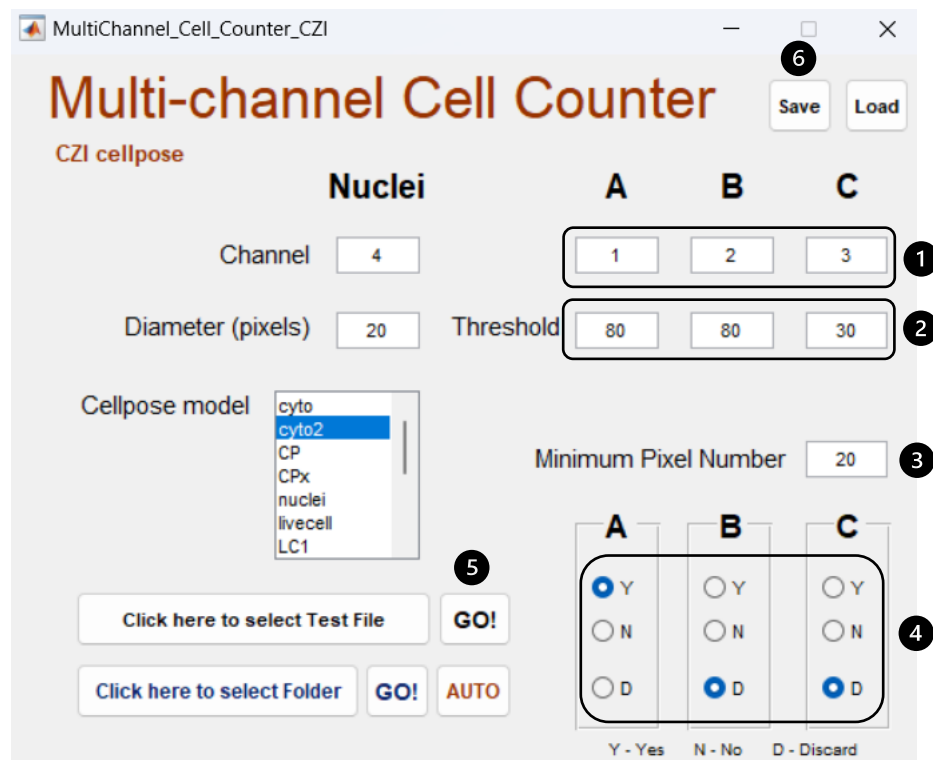

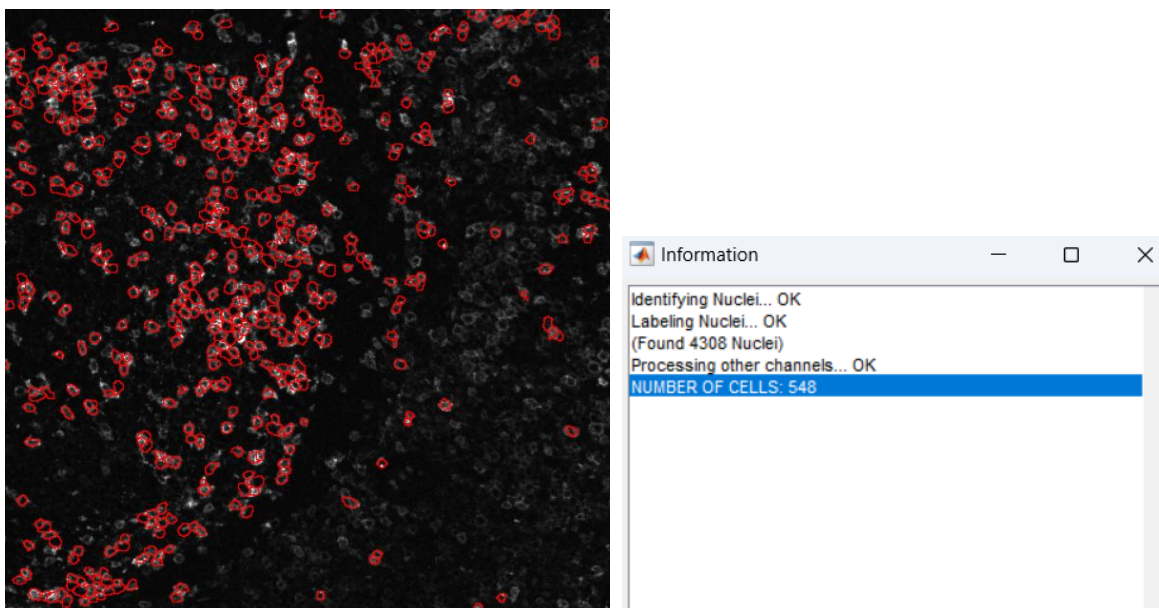


2 - Positive Cells for channel A - A(Y) B(D) C(D) - using Minimum Pixel Number 20

6. (Optional) Save processing parameters in an Excel file to be reused later.

Step 3 Process Folder

1. (Optional) Load processing parameters from Excel file
2. **Select Folder** with CZI files to be processed
3. Press **GO!**

CZI files in the folder will be processed using the combination set in A B C. An Excel file named MCC_results_[combination].xls and individual TIF files with positive nuclei outlines will be created.

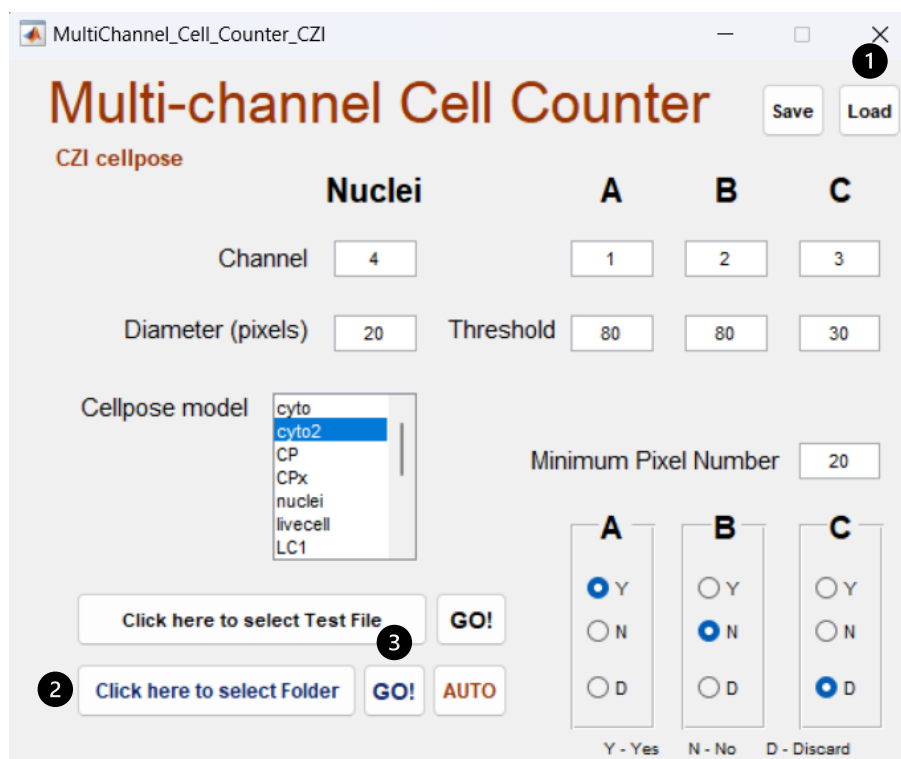

Step 4

Automatic Processing with different combinations

1. (Optional) Load processing parameters from Excel file
2. **Select Folder** with CZI files to be processed
3. Press **AUTO**
4. Select combinations to be processed and press **GO!**

CZI files in the folder will be processed using the combinations set in AUTO. Excel files named MCC_results_[combination].xls and individual TIF files with positive nuclei outlines will be created.

The image shows two windows from the software. The main window, titled 'MultiChannel_Cell_Counter_CZI', has a title bar with a close button (X) labeled '1'. The window contains the following elements: a 'Save' and 'Load' button pair; a 'Nuclei' section with a 'Channel' dropdown set to '4'; a 'Diameter (pixels)' input set to '20'; a 'Threshold' section with three inputs (80, 80, 30) for columns A, B, and C; a 'Cellpose model' dropdown menu with options: cyto, cyto2 (selected), CP, CPx, nuclei, livecell, and LC1; a 'Minimum Pixel Number' input set to '20'; a 'Click here to select Test File' button; a 'GO!' button labeled '3'; a 'Click here to select Folder' button labeled '2'; and a 'GO!' button labeled 'AUTO'. Below these are three columns of radio buttons for 'A', 'B', and 'C', each with options Y (Yes), N (No), and D (Discard). A legend at the bottom indicates 'Y - Yes', 'N - No', and 'D - Discard'. The second window, titled 'AUTO_sel...', has a title bar with standard window controls. It contains an 'ALL' button, a 'NONE' button, a grid of checkboxes for combinations of Y, N, and D for columns A, B, and C, and a 'GO!' button labeled '4'. The grid shows the following checked boxes: YYD, NYD, DND, NDY, and YDD.

NOTE

High DPI scaling issue

- If the graphical user interface (GUI) is not displayed as depicted in this Quick User Guide, you may need to override High DPI scaling in your Windows computer. To do so, right-click the **Multichannel_Cell_Counter_CZI** shortcut and select **Properties**. Click on the **Compatibility** tab and under Settings, select **Change high DPI settings**. In the High DPI scaling override section, select "Override high DPI scaling behavior. Scaling performed by:" and select **System (Enhanced)**.

HELP

For support, please contact joserino@medicina.ulisboa.pt